en **MARATHON®**
NO STING LIQUID SKIN PROTECTANT

Instructions for use

FOR EXTERNAL USE ONLY

Product Description:
Marathon® No Sting Liquid Skin Protectant is a non-cytotoxic (not harmful to the body) cyanoacrylate based, rapid drying, liquid barrier film for the protection of the skin. It is applied as a liquid and dries, within a minute, adhering to the contours of the skin to form a transparent flexible film. Marathon® No Sting Liquid Skin Protectant will wear off, naturally, as the skin regenerates.

Intended Use:

• Marathon® No Sting Liquid Skin Protectant is intended to protect intact or damaged skin from the effects of moisture, friction, (rubbing) and shear (tearing).

Indications For Use:

- Marathon® No Sting Liquid Skin Protectant helps protect skin exposed to irritation from moisture such as urine, feces, digestive juices, perspiration and wound drainage.
- Marathon® No Sting Liquid Skin Protectant can also be used in areas that are exposed to friction and shear from bedding, clothing, shoes, or any material that would rub against the skin.
- Marathon® No Sting Liquid Skin Protectant helps protect the skin against irritation caused by adhesive products.


Contraindications:

- Do NOT apply directly to deep, open, or bleeding wounds.
- Do NOT apply to chronic wounds (persistent or unhealing wounds).
- Do NOT apply to second or third degree burns.
- Do NOT apply to infected areas.
- Do NOT use on individuals with a known sensitivity to cyanoacrylates or formaldehyde.

Instructions For Use:

- Clean and dry the skin (Remedy Antimicrobial or Foaming Cleanser is suggested) before applying Marathon® No Sting Liquid Skin Protectant. Make sure that the skin is thoroughly dry and that all traces of moisturizers, lotions, lotion soaps or lanolin are removed from the application site.
- Protect any sensitive areas, such as stomas and mucous membranes (e.g. the gums and inside of mouth).
- Remove the applicator from the packaging. While holding the applicator upright (with the sponge tip to the top), firmly squeeze between the finger and thumb to crush the sealed inner tube.

Correct Activation:

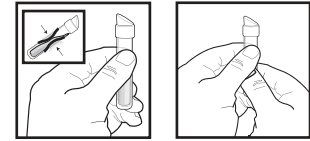One handed method

When two hands are needed

INCORRECT ACTIVATION:

DO NOT try to bend or snap the tube to activate.

- Turn the applicator upside down and wait for 5 seconds while the Marathon® No Sting Liquid Skin Protectant soaks through the sponge tip. Gently squeezing the tube will speed up this process, but be cautious of over squeezing.
- As soon as the sponge is wet, slowly spread the Marathon® No Sting Liquid Skin Protectant, gently and evenly, to form a thin transparent coating. Ensure that the desired area is completely covered and extends at least an inch (2,5cm) beyond the affected region.

- If Marathon® No Sting Liquid Skin Protectant is used at room temperature it will dry within a minute. Some people may feel a temporary warming sensation as the cyanoacrylate dries.
- Do not over apply or rub the Marathon® No Sting Liquid Skin Protectant in like a cream as the coverage will not be as effective and the warming sensation will be increased. It is most effective when the coating is paper-thin.
- Ensure that Marathon® No Sting Liquid Skin Protectant has dried, completely, before covering with clothes,

applying adhesive products such as tapes and dressings, or incontinence products.

- Marathon® No Sting Liquid Skin Protectant naturally wears off over 1-3 days depending on skin type and activity level. Inspect the film regularly (at least once a day) and if any cracking or flaking appears, apply another layer on top, until the area is completely covered once again. It is not necessary to remove the cracked and flaking layer, however, if it is applied too frequently excess build-up will lead to a thick inflexible film.
- One Marathon® No Sting Liquid Skin Protectant applicator will cover an area of approximately 4 in x 4 in (10 cm x 10 cm).
- Do not use cleansers that leave a protective barrier, such as dimethicone, petrolatum, mineral oil, etc.

Warnings:

- If irritation (redness, inflammation, heat, pain, and change in skin color) develops, immediately discontinue product use. If pain, and change in skin color develops, immediately discontinue product use. If irritation persists, seek medical attention.
- Only use for up to 30 days on damaged skin.
- For external use only.
- Do not reuse in whole or part as it may compromise the performance of the protectant.
- If Marathon® No Sting Liquid Skin Protectant is applied to a moist skin surface it will set very rapidly and a noticeable, but temporary warming sensation will be felt.
- Do not use under drug delivery patches.
- Avoid contact with mucous membranes such as the eyes and mouth.
- Keep out of reach of children.

Precautions:

- Avoid contact with materials such as latex gloves, instruments, clothing, bedding until cyanoacrylate has set (approximately 1 minute).

Adverse Reactions:

- If Marathon® No Sting Liquid Skin Protectant enters the eye immediately flush with lukewarm water. Do not rub, or try to separate the eyelids. Seek medical attention.
- If Marathon® No Sting Liquid Skin Protectant accidentally contacts an open wound, flush with warm tap water. Removal is not necessary.
- If the skin or other materials become stuck to the area, gently roll apart to separate. If unable to separate easily, seek medical attention.
- If pain occurs, and persists, seek medical attention.

Removal:

- Excess Marathon® No Sting Liquid Skin Protectant can be removed before it dries by wiping with a clean tissue. Once it has dried, unwanted cyanoacrylate can be softened with petroleum jelly and removed. It will naturally wear off within 1-3 days as the skin renews itself.

Storage:


- For best results, store Marathon® No Sting Liquid Skin Protectant between 5°C and 25°C in a cool, dry place. Protect from freezing and avoid excessive heat and humidity.
- For shelf life refer to expiration date on product.

Contents:

- Each applicator contains 0.5g of a cyanoacrylate blend formulation.
- Each box contains 5 applicators of Marathon® No Sting Liquid Skin Protectant.

Disposal:

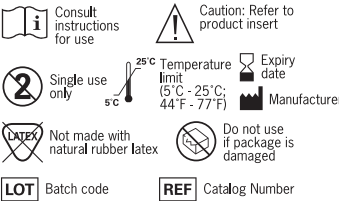- Throw away the applicator into standard waste after use.

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Medline and Marathon are registered trademarks of Medline Industries, Inc. - www.medline.com/woundcare
Medline Industries, Inc., Mundelein, IL 60060 USA
Distributed in Australia by: Medline International Two Australia Pty. Ltd., Sydney, Australia